


Unistik® 3

Pre-set single use safety lancets with *Comfort Zone Technology*®

When pressed against the sample site, a series of 8 raised dots on the platform stimulate the nerve endings. This sends a powerful message of comfort to the brain – masking the weaker pain stimulus from the lancet.

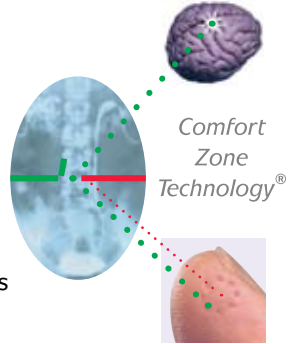

The Unistik 3 range of single use safety lancets has 4 models for use across all clinical environments.

Unistik® 3 Comfort

Gauge: 28G. Depth: 1.8mm
Average Blood Volume: 5µL-30µL.
Suggested Test Areas e.g.
Ward blood glucose monitoring. HbA1c.
Suggested Patient Areas e.g.
Diabetes Clinics. Paediatrics. Geriatrics.
Home Use.



Box of 100: AT1042
Box of 200: AT1044

Unistik® 3 Normal

Gauge: 23G. Depth: 1.8mm
Average Blood Volume: 30µL-75µL.
Suggested Test Areas e.g.
Ward blood glucose monitoring.
Cholesterol. HbA1c.
Suggested Patient Areas e.g.
Adults. Geriatrics. Home Use.
Blood transfusion centres.



Box of 100: AT1002
Box of 200: AT1004

Unistik® 3 Extra

Gauge: 21G. Depth: 2.0mm
Average Blood Volume: 75µL-125µL.
Suggested Test Areas e.g.
Cholesterol. Haemoglobin (Hgb). Diagnostic kit tests such as cardiac enzymes.
Suggested Patient Areas e.g.
Adults (tougher skin). Blood transfusion centres.



Box of 100: AT1012
Box of 200: AT1014

Unistik® 3 Neonatal & Laboratory

Gauge: 18G. Depth: 1.8mm
Average Blood Volume: 200µL+.
Neonatal Suggested Test Areas e.g.
Congenital Hypothyroidism. Phenylketonuria (PKU). Single heelstick tests such as Bilirubin, Glucose, Haemoglobin. Sickle Cell disease. Other diseases.
Neonatal Suggested Patient Areas e.g.
Heelstick samples on neonates.



Box of 100: AT1052
Box of 200: AT1054

Laboratory Suggested Test Areas e.g.
Blood Gases & Electrolytes. Chemistry Test Profiles. H.pylori.
Laboratory Suggested Patient Areas e.g.
ITU / ICU. CCU. Laboratories

How to use Unistik® 3

...simply twist


Hold the Unistik 3 body and twist off the **lancet cap** until you feel it separate from the device.



Don't pull, just twist. Dispose of the cap in a suitable container.

...then sample

Choose a site as indicated by the shaded areas on pictures **A** and **B**. Press the *Comfort Zone Technology* **platform** firmly against the chosen site and press the **release button**.

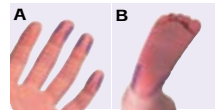

If required, massage the sample site, taking care not to squeeze too hard at the site.



The single use needle retracts immediately after sampling, leaving the device safe for immediate disposal. Use a suitable sharps container.

 **OWEN MUMFORD**
www.owenmumford.com

AT/1042/MCC/1205/8

UK / EU

Owen Mumford Ltd
Brook Hill, Woodstock
Oxford OX20 1TU, United Kingdom
T : 01993 812021
F : 01993 813466
E : info@owenmumford.net

NORTH AMERICA

Owen Mumford Inc
1755 West Oak Commons Court
Marietta, Georgia 30062, USA
T : 800-421-6936
F : 770-977-2866
E : info@owenmumfordinc.com

FRANCE

Owen Mumford
BP 444
27204 Vernon Cedex, France
T : 02 32 51 88 70
F : 02 32 21 96 95
E : contact@owenmumford-fr.com

GERMANY

Owen Mumford GmbH
Alte Häge 1
63762 Großostheim-Ringheim, Germany
T : 06026/97750
F : 06026/977529
E : info@owen-mumford.de

EXPORT

Owen Mumford Ltd
Brook Hill, Woodstock
Oxford OX20 1TU, United Kingdom
T : +44 1993 812021
F : +44 1993 813473
E : export@owenmumford.co.uk